

Eclipse[®]

Superabsorbent Dressings



Exudate management: critical in wound healing

Well managed exudate can¹:

- Shorten wound healing time
- Optimise resources
- Reduce the impact on the patient, both physically and psychologically

Eclipse superabsorbers have Advancis' CrystaLock technology:



- Handle high levels exudate
- Lock in core
- Reduce leaks
- Reduce risk of maceration


Dynamic fluid uptake maintains moist wound healing

Dynamic removal of excess exudate:

- Ensures a moist wound healing environment is established and maintained


Optimal healing environment with the removal of detrimental cells like:

- Infection causing bacteria
- MMPs that can stall wound healing

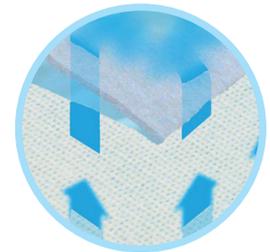

Zero strikethrough to improve patient Quality of Life

- Frequent dressing changes and exudate leaking onto skin or clothes, can be distressing for patients, families and carers
- Patients may experience social isolation and anxiety due to leakage of dressings²
- Eclipse is the only superabsorber with zero strikethrough³

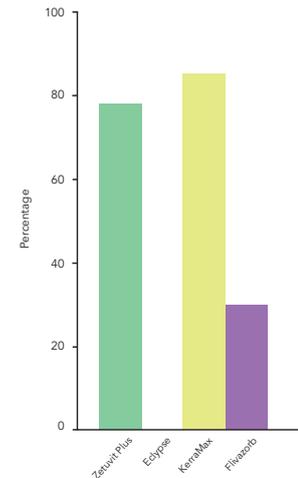

Strikethrough at dressing change on day 7

Unparalleled range and unique shapes

The **only** range of superabsorbers to feature standard, soft adherent with Silfix technology and bordered versions.

Plus unique shaped dressings are available with Eclipse Foot, Eclipse Boot and Eclipse Contour.



NEW
10x20cm and
20x20cm sizes

As the largest range available, you can be confident that there is an Eclipse for all you patients!

Four layers of confidence for effective exudate management



Rapid wicking layer

- Draws exudate from the wound bed
- Ensures the wound doesn't dry but maceration is avoided
- Exudate is delivered into the **CrystaLock** layer



Sub-backing layer

- Helps retain the dressing's shape, encouraging even distribution of exudate



Backing layer

- Non-strikethrough, preventing the patient embarrassment associated with leaks
- High moisture vapour transfer rate (MVTR)
- Water-resistant
- Bacteria- and viral-proof



CrystaLock layer

- Absorbs exudate and locks it within the dressing to prevent leakage
- Associated bacteria is locked away too, which reduces the risk of infection



Exudate: Managed

Venous leg ulcer

Ulcer size had increased due to poor exudate management

Eclipse[®] Contour used under compression bandaging


Significant improvement at week 4 with some areas healed



Diabetic foot ulcer

Following debridement, **Eclipse[®]** used to manage high exudate levels and prevent maceration of intact tissue

Wounds healed or almost healed within 8 weeks

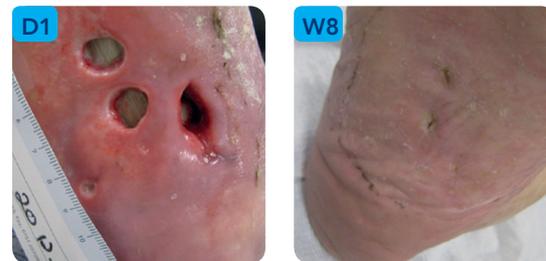

Eclipse® Superabsorbent Dressings



Ordering information

Eclipse®



Size	Stock code	Pack size	NHS code	PIP code
10x10cm	CR3818	20	EJE029	372-1057
15x15cm	CR3769	20	EJE032	305-1349
10x20cm	CR4460	10	coming soon	coming soon
20x20cm	CR4439	10	coming soon	coming soon
20x30cm	CR3743	20	EJE030	306-1272
60x40cm	CR3808	10	EJE034	232-1131

Eclipse® ADHERENT



Size	Stock code	Pack size	NHS code	PIP code
10x10cm	CR3881	10	EKH034	325-2293
10x20cm	CR3883	10	EKH011	325-2269
15x15cm	CR3863	10	EKH036	325-2277
20x30cm	CR3864	10	EKH015	325-2285

Eclipse® ADHERENTSACRAL



Size	Stock code	Pack size	NHS code	PIP code
17x19cm	CR3985	10	EME043	348-2213
22x23cm	CR3986	10	EME044	348-2239

Eclipse® BORDER



	Size	Stock code	Pack size	NHS code	PIP code
OVAL	15x15cm	CR4296	10	EKH061	383-6897
	20x30cm	CR4297	10	EKH062	383-6905
	10x20cm	CR4300	10	EKH063	383-5881
	15x20cm	CR4299	10	EKH064	383-5899
	20x30cm	CR4298	10	EKH065	383-5907

Eclipse® CONTOUR



Size	Stock code	Pack size	NHS code	PIP code
30x51cm	CR4398	10	EJE166	400-6219

Eclipse® FOOT



Size	Stock code	Pack size	NHS code	PIP code
33x48cm	CR4219	5	EJE086	376-5468

Eclipse® BOOT



Size	Stock code	Pack size	NHS code	PIP code
55x47cm	CR4237	5	EJE085	380-4630
60x70cm	CR4014	5	EME042	347-6207
71x80cm	CR4229	5	coming soon	408-0750